


MANUAL

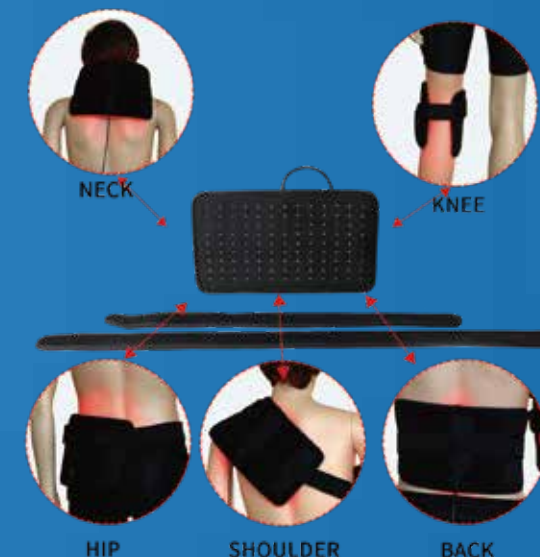

Please read the manual carefully before using the device

Attention

- ⚠ People with heart disease or heart pacemaker can not use this light therapy system.
- ⚠ Pregnant woman, skin-ulcers and allergic people can not use it.

Features

Model:	RLT-WB01
Material:	75% neoprene and 25% velvet
Effective area:	1pcs 12 X 6inch
Power supply:	AC100-240V or DC 12V DGYAO Portable Battery
Frequency:	50/60 Hz
Wavelength:	660nm+880nm
Number:	52pcs Red lights & 65pcs IR lights
Light Beads Craft:	DIP
Temperature :	Max 104 °F
N.W:	530g
G.W:	1120g

Contact Us

- Customer Service Number: +86-139-2922-1675(skype)
- Email: dgyao@dgyao.net
- Customer Care Hours: 6am to 11am, 6pm-12pm
- Standard Time Monday – Sunday Except Major Holidays

DONG GUAN TOP YAO INDUSTRY CO.,LTD

- ▶ Email: dgyao@dgyao.net
- ▶ Facebook: [@dgyao](https://www.facebook.com/dgyao)
- ▶ Skype: +86-139-2922-1675
- ▶ Amazon Store: <https://www.amazon.com/dgyao>
- ▶ Official Website: <https://www.dgyao.net>



www.dgyao.net



FACEBOOK



TWITTER



YOUTUBE



LINKEDIN

Warranty and disclaimer

Limited One-Year Warranty

DGYAO carries a limited one-year warranty for the original purchaser against any technical defects. Warranty coverage applies to the original owner and to the original product only and is not transferable.

DGYAO warrants that the device shall be free from defects in material and workmanship for one year from the date of purchase when the device is used and maintained according to the requirements outlined in the User Manual.

If your device stops functioning or working as intended during the one-year warranty period, please call DGYAO Customer Care at Skype +86-139-2922-1675 and one of our customer service team members will work with you to understand the issue and give you instructions on how to get your DGYAO LIGHT THERAPY repaired or replaced, free of charge.

This warranty does not cover normal wear of the device and does not apply to any device that has been tampered with or damaged due to misuse, abuse, or negligent handling.

Dear Customer, thanks so much for your great support on us! If there is anything problem with device, please feel free to contact us and we will be glad to help you.

Please send us the photos or videos of the product failure, and please kindly advise me your order number, phone number, contact name and shipping address.

If your order is still on one-year warranty, we will resend you the new replacement once we confirm the exact problems.

Thanks for your cooperation in advance.

Specific Warranty Policy:

- ✓ 1-Year defects liability period
- ✓ 60 days return and replacement warranty
- ✓ 18-month free replacement for charger, controller and connectors. (if they are broken).
- ✓ Return postage borne by customers
- ✓ Return requests as “I don’t like it”, “I don’t want it”, “the color doesn’t look good” are not accepted.
- ✓ The-less-than-three-light-broke-situation will not covered in the warranty policy because it doesn’t affect usage effect.
- ✓ keep the original outer package for successful return.

Frequently Asked Questions

Q1: What's the wavelenath of the lights?

A1: 660nm in red lights and 880nm in infrared lights.

Q2: Does it has both red and infrared lights?

A2: With 65pcs near infrared lights + 52pcs red lights, total 117pcs.

Q3: Can I use this device in Europe?

A3: Yes, you only need a EU adaptor and it can be used in Europe.

Q4: Does it need to be plugged into the socket?

A4: Yes. There is also another model that is directly connected to the DGYAO power bank and can be used similarly(Sold separately).

Q5: Can multiple people share the device?

A5: Yes.The device can be used by different people and your pets too if you like.

Q6: What voltage can this device be used?

A6: The device can be used within AC110-240V.

Q7: Is it more effective on bare skin than over cloth?

A7: Yes, it is more effective on bare skin than over cloth.

Q8: Are there side effects?

A8: Don't worry about it, but please notice you may will have skin irritation if you 're allergic to the fibers.

Q9: How long will it take to experience relief of device?

A9: We recommend wearing it for 2-3 therapy sessions perday for the first month.20 minutes each time.Take arest and drink some water to relax yourself before next following use.Please be patient because it may take several weeks to experience optimal relief.

A10: Can people with heart disease or heart pacemaker use it?

A10: We recommend that you consult your primary care provider or physician when introducing Light Therapy.

- Contact your Doctor if you have any medical questions regarding use of this device.Contact your Doctor prior to using-taking any light.sensitive medications or medical questions.

- Light Sensitivity Test- the User Guide also recommends testing yourself for light sensitivity priorto initial treatment.

Q11: What situations should I not use the therapy?

A11: Pain which may benefit from heat can be relieved by red and infrared light therapy.But under some situation,do'NOT use the light therapy:

- 1.High Fever.Infrared light may increase the temperature of your tissue and can worsen fever.
- 2.Tumour &tuberculosis.Improper use of infraredradia tion may lead to more severe tumors.
- 3.Patients fitted with heart pacemaker can NOT use light therapy.
- 4.Acute soft Tissue Injury. Bruises, sprains,bruises, falls or sport injury cause subcutaneous soft tissue torn bleeding or exudation. Do Not use light therapy within 24-48 hours, otherwise it will increase inflammation and swell.

For more information, please go to the blog of our official website: www.dgyaodevice.com

Safety Requirements and Notes

Dear Users, thanks for purchasing red light & Infrared light therapy. Instruction for use has been provided to ensure that you will be able to understand and safely use this device.

This instruction illustrate the device and shows you how to operate and care for the unit. Please read thoroughly before using your new unit, and follow the instructions.

Important Notice


If you have any problem or need help in using, please contact our technical service center for help in time. We will give you technical support or arrange professional technical expert for service at the first time.

Use correctly can extend the life of equipment, also make fully usability of equipment by farthest.

Abnormal operation may do harm to equipment or personal safety. Our company is irresponsible for abnormal condition, hazards to equipment or personal injury caused by operations that are absolutely prohibited as specified in this document. Disclaims any responsibility for safety, reliability or performance of this equipment by not observing the instructions!

Any faults arising from such non-observance will invalidate the warranty!

Please carefully read the instructions in the "safety requirements", "Note" and the special warnings "!" part of the content.

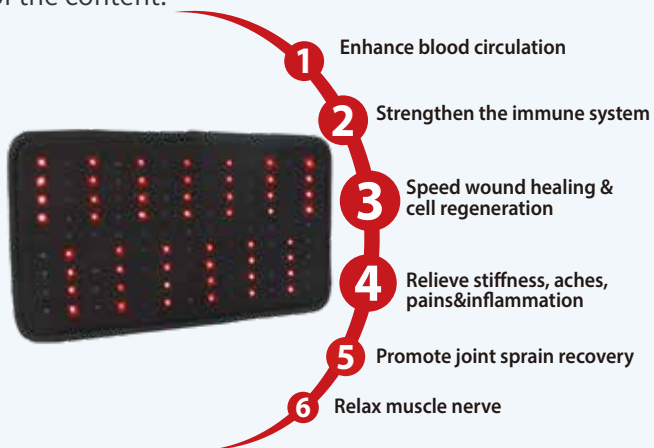


Introduction

Effective wavelengths of near infrared light therapy include 750 nm to 1,400 nm. The wavelength of the near infrared is more beneficial in treating deep wounds and relieving muscle and joint pain. 850nm and 980 nm are common for near infrared light therapy.

DGYAO adopts 65pcs 880 nm infrared led lights (invisible to naked light) and 52pcs 660 nm red led lights. The combination of red and infrared light increase blood circulation by thermal effect to provide pain relief.

How does Infrared light work?

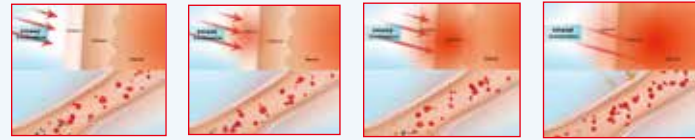

The NIR is the shortest wavelength and the most penetrating infrared ray, which can penetrate 5 to 10 mm of the human body tissues. Infrared light can penetrate the skin, directly to the muscles, subcutaneous tissue which helps to accelerate blood circulation, increase metabolism, relieve pain and increase muscle relaxation, resulting in massage effects. Infrared treats the disease at different levels because of mobilize the body' s own disease-resistant ability.

Function

SAFE & NATURAL PAIN TREATMENT:

- ★ Relieve stiffness & Spasm
- ★ Reduce pain & inflammation
- ★ Safe for all types of skin
- ★ Enhance recovery from exercise or injury
- ★ Improve hydration & oxygenation of the tissues
- ★ Relax muscle cramps & bone joint Torn Meniscus pain
- ★ OTC class II medical device and LED lights reach deep into joints and soft tissue
- ★ Improve skin texture
- ★ Calm the nervous system
- ★ Easily administer at home

DGYAO Infrared & Red Light Therapy Device penetrates muscles and joints effectively to improve blood circulation and relieve pain.




Product description

LED red light consists of:

- ① 1 X Infrared light Pad
- ② 117pcs Lights: DIP, 9W (52pcs Red lights & 65pcs IR lights)
- ③ 1 X Auto shut-off timer: 20 Minutes
- ④ 1 X AC Adapter: AC 100~240V with a cord
- ⑤ 1 X DC Cord


- Infrared light
- Red light

Red light: 52PCS
Infrared light: 65PCS
Total light: 117PCS

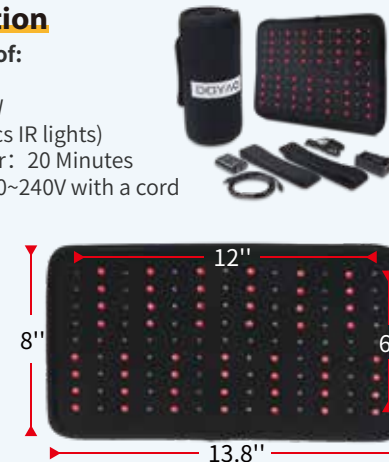

Instructions for use

- ★ A1: 1 X Red & Infrared LED Lights Therapy Devices
- ★ A2: 1 X DC Cord.
- ★ A3: Auto Shut-Off Timer (20 Minutes).
- ★ A4: 1 X AC Adapter (AC 100-240V) & Power Cord.

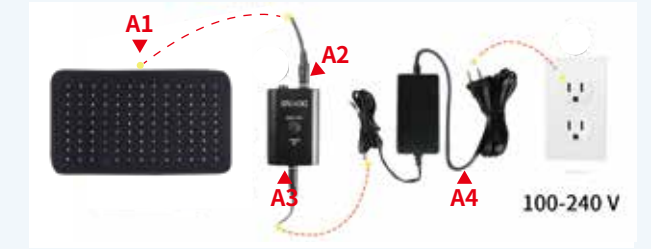

AUTO SHUT-OFF TIMER (20 MINUTES). Do NOT exceed 40 minutes of use in one therapy session.

DC cord and timer connection



The DC cord connects to the plug **OUT1** or **OUT2** of Timer.



The AC power adapter connects to the plug **IN** of Timer.

Cleaning

1. Disconnect the device from the power supply before cleaning.
2. Using a soft cloth with a little water to gently wipe the surface of device.

Caution

1. Do not use alcohol or other corrosive cleaning agents to avoid corroding the surface and damaging internal parts.
2. Do not apply liquid cleaner directly to any surface of the adapter, timer, pad as it could cause damage to the electrical components.

Storage

The device must be stored in dry area; away from high humidity;

Warning

- No washing
- No compression
- Stay dry
- Handle lightly