



USER MANUAL

Cordless Joint Wrap



RLT-MS23-KEP

Email: dg Yao@dgyao.net



01 | PRODUCT STRUCTURE




- 1 Knee Pad
- 2 Magic Glue
- 3 LED Lights
- 4 LOGO
- 5 Magnetic interface
- 6 Controller Battery
- 7 Extension Strap
- 8 Storage Bag
- 9 Type-C Cable
- 10 Power Plug

02 | OPERATION GUIDE



Set first! Then match power to the wrap, enjoy it! **There are 5 magnetic points, they must be set to the same points on belt.** Make sure a clear tick sound on each magnets are sound.

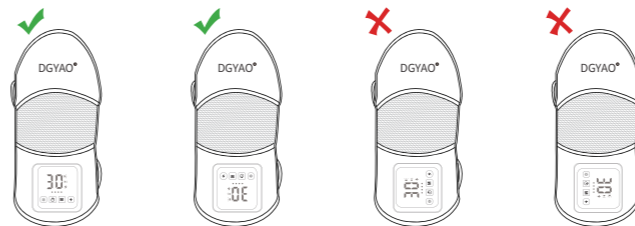

03 | CONTROLLER BATTERY

If the controller has not been used for a long time, please charge the controller in about 2 months!



- Three Lighting Mode
- Five Light Intensity
- Pulsing Mode
- 5-30 Mins set

- Power:** Long press 3 seconds Power On/Power Off.
- Intensity:** Press to change light intensities. H1: 40% H2: 55% H3: 70% H4: 85% H5: 100% full intensity. The higher the intensity, the brighter the light.
- Pulse:** Click the strobe to turn on, the icon is on, and the frequency is fixed at 10HZ.
- Modes:** Press "M" to set modes levels
 ◀▶: Red lights+Infrared lights ◀: Red lights only ▶: Infrared lights only
- Timer Set:** Press to set auto off timer set at a range from 5-30mins. 5 mins per press. When the countdown ends, the device will beep 3 times, and then the device will automatically shut down.

Screen will show Real-Timing 10 seconds after intensity set. Display will auto off to enter power saving mode after 3 minutes no action. Press power button will wakeup the screen again.

04 | EASY-START SET GUIDE

GET STARTED

PREPARING TO USE

- Please charge the power supply before your first use.

1

- Clean the magnets and the lights. Check all the condition of whole unit.

POWER ON AND CHOOSE SETTINGS

- Long press “☺” 3 seconds to power on the controller.
- Choose your settings of timer, intensities, pulse rate.

2

CONNECT BATTERY TO BELT

- Set the power supply to the belt.
- Each magnets must be matched in the same points.
- Wrap the whole belt on your body and enjoy it now!

3

HOW TO CHARGE:

Indicator Lights 

1. Turn off the battery. Use the cable and charger plug to charge the battery.
2. White indicators will show actual battery level.
3. Indicators will slowly flash to show charging status.
4. Controller will auto lock when charging. So you can not use it when charging.
5. It takes about 3 hours to fully charged the battery.
6. When the screen shows "LO" with a beep, it means that the battery is low and needs to be charged.

05 | SPECIFICATIONS

Model:	RLT-MS23-KEP
Weight:	1.1LBS / 0.5KG
Dimensions:	12.2" × 7.2"
Battery Capacity:	4000 mAh
Power Density:	Max 80 mW/cm ² @ 1cm
Low EMF:	<0.1uT @ 0cm
Wavelength:	660nm + 850nm
Pulse Rate:	10HZ
Timer Set:	5-30 minutes
Input:	DC 5V/3A
Lights:	96 Red lights+ 48 infrared lights
Light Beads Craft:	SMD 3 IN 1 Bulb
Temperature:	Max 104°F
Material:	75% neoprene and 25% velvet

06 | WARRANTY

THANK YOU FOR CHOOSING DGYAO

The dedication we put into engineering high-performance therapy devices extends into after-sales service. If you face any problems with your DGYAO product, DO NOT return it to your retailer. Get immediate assistance from DGYAO Customer Support via email to dgyao@dgyao.net.

- 1-Year defects liability period.
- 60 days return and replacement warranty.

If DGYAO THERAPY ITEM ever breaks, we'll replace it—free of charge.

05 | TROUBLE SHOOTING GUIDE

Quickly troubleshoot simple issues you might experience with Portable Therapy Wrap using the table below.

Please contact a customer support representative for additional assistance at: dgyao@dgyao.net.

ISSUE	POSSIBLE REASON	HOW TO RESOLVE
DEVICE IS NOT WORKING	Magnets are matching in opposite way.	Take off the battery and check the match. Each magnet must be set in certain direction.
	Device is not turned on.	Press the central button to turn on the device.
DEVICE STOPPED WORKING SUDDENLY	Auto off time is up. It auto shut off.	Turn the system on to begin a new red light therapy session.
LIGHTS TURN OFF WHEN USING	Battery is out. Needs to be charged.	Please charge the battery with our charger.

07 | CARE INSTRUCTION

1. Do not get the device near water or fire.
2. This belt is not machine washable or dryer usable.
3. Simply clean it by wipes or towels which is not very wet.
4. The lights will function still if having some sweat on your back.

07 | CONTACT US

DONG GUAN TOP YAO INDUSTRY CO.,LTD

- Email: dgyao@dgyao.net
- Facebook: [@dgyaocn](https://www.facebook.com/dgyaocn)
- Skype: +86 139 2922 1675
- Amazon Store: <https://www.amazon.com/dgyao>
- Official Website: <https://www.dgyao.net>

Consumers and clinicians call: +86-139-2922-1675 (Skype)

Customer Care Hours: 6am to 11am, 6pm-12pm Eastern Standard Time Monday - Sunday Except major holidays.



@dgyaocn
FACEBOOK



@DGYAO
YOUTUBE



+13929221675
WHATSAPP



@dgyaocn
INSTAGRAM



@dgyao_rlt
TIKTOK